


## **Table Top Autoclave Class B Series**







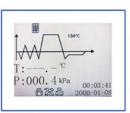







## Features:

\* Built-in open type water tank

The sterilizer adopts easy-clean open type water tank which can support repeated program running if fully injected with water.

\* High-efficiency ultimate vacuum

The sterilizer adopts high-efficiency low noise vacuum system which has excellent effects.

\* Large LCD display for BKMZA/BKMZB

The LCD screen can display temperature, pressure, time, operating status, failure warning and other information. It is convenient for customers to observe the sterilizer running status.

\* Multiple program types

The system has various programs that include packed items, unpacked items, BD testing program, vacuum testing program and drying function.

Custom program, rapid program and preheat function (for BKM-Z16B).

\* Standard USB port for BKMZA/BKMZB

Users can store the sterilization data with USB disk.

\* Optional mini printer can be attached to record the process of sterilization.

## **Technical Parameters:**

Model	BKMZA (Previous name: BKM-Z18B)	BKMZB (Previous name: BKM-Z24B)	
Capacity	18L	23L	
Chamber Size(mm)	φ247*360	φ247*470	
Sterilization Class	Class B		
Sterilization Temperature	121℃, 134℃		
Drying System	Vacuum drying system		
Display	LCD display		
Testing System	B&D Test		
	Vacuum Test		
	Helix Test		
Control Precision	Temperature: 1℃		
	Pressure: 0.1bar		
Sterilization Data	BKM-Z16B: Printer(optional)		
	BKMZA/BKMZB: USB(standard) and printer(optional)		
Safety System	Hand lock door		
	Pressure lock system		
	Relief valve in case of over pressure		
	Pressure or temperature over load protection		
	Alarm for system failure, finish reminding, water level warning		
Water Supply System	Build-in water tank easy to clean		
Water Tank Capacity	3.5L		
Water Consumption	0.16L~0.18L in one cycle		
Tray Holder	3 pcs SS trays on SS shelf		
Chamber	SUS304		
	Max working pressure: 2.3bar		
	Min working pressure: -0.9bar	Min working pressure: -0.9bar	
	Design temperature: 140℃		
Ambient Temp.	5~40℃		
Noise	<65dB		
Consumption	1950VA	1950VA	
Power Supply	Standard: 220V±10%,50Hz/60Hz; Optional: 110V±10%,50Hz/60Hz. Wiring Method: Plug		
External Size(W*D*H)mm	495*600*410	495*700*410	
Package Size(W*D*H)mm	610*810*590	610*810*590	
Gross Weight(kg)	63	65	

Note: It is applicable to the altitude below 3500 meters, and the parts need to be changed when the altitude is over 2700 meters.

03 04